



Adjustable
Trolli

ORMESA®
• MADE IN ITALY SINCE 1980 •

More than
a folding
system!



Practicality
at hand



Advanced Design for a no limits mobility

Adjustable Trolli by Ormesa is a compact and light **seating system for teens and adults**, ideal to move in any environment.

It is easy to handle thanks to the **international patented system** which allows you to **fold it with one single operation**. Once folded, it remains compact, **occupying the minimum space**.

Adjustable Trolli has a **stable and resistant structure** guaranteed by the **double cross** positioned behind the backrest and under the seat which gives greater stability and comfort.

The **upholstery** is **made with advanced technological materials**: it is extremely **breathable, resistant, easily washable** and **fire retardant**. It is available in **two colors**: red and black.

Thanks to a **wide range of available components**, it easily meets every need.



Adjustable Trolli is CRASH TESTED according to **ISO 7176/19** and **ANSI RESNA WC/19** standards.

PRODUCT
MARKED **CE**

ADJUSTABLE TROLLI is a MEDICAL DEVICE COMPLYING WITH THE SAFETY REQUIREMENTS DEFINED BY REGULATION (EU) 2017/745



Practical and light

Sea, mountain, city: no problem!




Adjustable Trolli has an aluminium **folding frame**, painted with non-toxic epoxy powders.

Once folded it remains compact: this allows you to **carry it on its rear wheels like a travel trolley**.

Its **low weight** and its **simplicity of folding**, makes it easy to lift or move: it is perfect if you need to store it in the car or if you need to take public transport.

Adjustable Trolli is the **ideal solution** for those looking for **stability, practicality and comfort, all in one!**

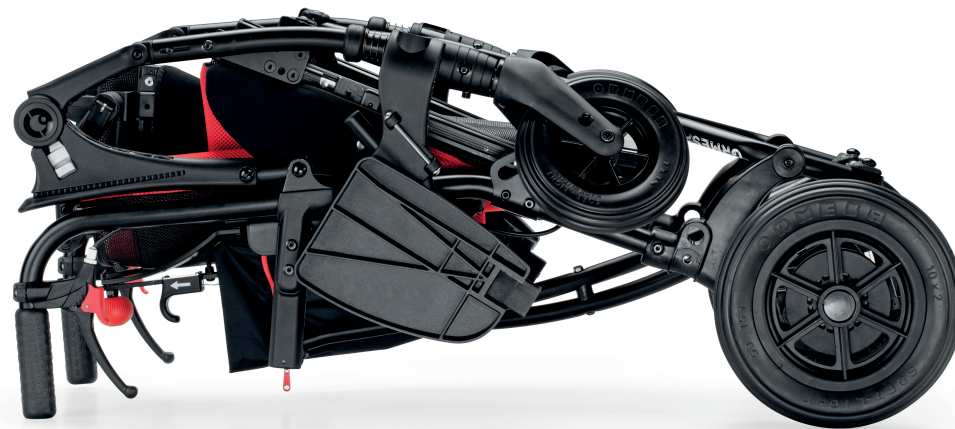

Technical features and adjustments

Upholstery

The upholstery of Adjustable Trolli has a **honeycomb** leading-edge structure, therefore it is **extremely breathable in contact with the skin**. The structure of the fabric creates an “air chamber” that allows both **air circulation for maximum breathability** and an **effective natural thermal insulation**.

The double layer is a guarantee of:

- › **Hygiene:** Thanks to the excellent breathability, it reduces sweating and bacteria formation.
- › **Washing:** despite the upholstery is not removable, the fabric composition allows easy sanitization and washing.
- › **Quick drying.**
- › **Resistance:** High resistance to abrasion. The high resistance texture maintains its performance even at high loads (ISO 12947 test).
- › **Fire retardant:** according to EN 1021-1:2014 and EN 1021-2:2014.





Postural headrest

Adjustable in tilt and depth thanks to levers and height adjustable by screw. Safe adjustment without the risk of losing the screw. It closed with the stroller, there is no need to remove it.



Backrest

Reclining in three positions. Simple adjustment in height by sliding, without the risk of losing the screw. With a bag on the back to store the manual and the supplied wrench. The back cross ensures stability, comfort and endurance.



Seat

Ergonomic and slightly tilted seat to prevent slipping of the pelvis and ensure greater comfort. Adjustable in depth by sliding and in width thanks to the additional component 958 padded pelvic side supports. Simple adjustment without the risk of losing the screw. The cross below ensures stability, comfort and endurance.



Handles

Push handles with manual and parking brakes. Very useful handles to manage the stroller in different environments. Adjustment in height with the component 951.



Legrest

Adjustable in height.

Footrests

Tilting and swivel footrest in line with the frame without encumbrance. Adjustment in tilt with the component 956.

The calf support band is also included for more comfort.



Wheels

With "Ormesa Design" original profile. Solid, light and **all-terrain wheels**.

Front solid swiveling wheels. 7" (17 cm) diameter for size 36 and 40; 8" (20 cm) diameter for size 43.

Rear solid wheels with independent drum brakes, 10" (25 cm) diameter for size 36 and 40; 12" (30 cm) diameter for size 43. **Overcoming step bar** on both rear wheels.

Components



838 Trunk side supports
adjustable in height, width and rotation.
Available in two sizes.



947 Pelvic belt with variable angle
Available in two sizes.



920 Four point pelvic belt
Available in two sizes.



853 Vest harness with zipper



853 slim Four point shaped harness



834 Padded abduction block
removable. Available in two sizes.

Components



839 Front handlebar

Adjustable in tilt. Easy to remove to fold the stroller.



960 Heel rests



956 Footrests adjustable in tilt



858 Shopping basket



812 Front wheel directional locks



844 Padded side protection covers

Components



891 Set of 4 tie down hooks

for the forward-facing transport of the pushchair in a motor vehicle (private/cars, buses..).

See User and Maintenance Manual provided with the aid.




951 Height adjustable push handles

It increases the total height of max. 7 cm/2,7 in (94-101 cm/37-39,7 in)




Adjustable Trolli sizes and weights

Unit of Measurement: cm (in); kg (lb)

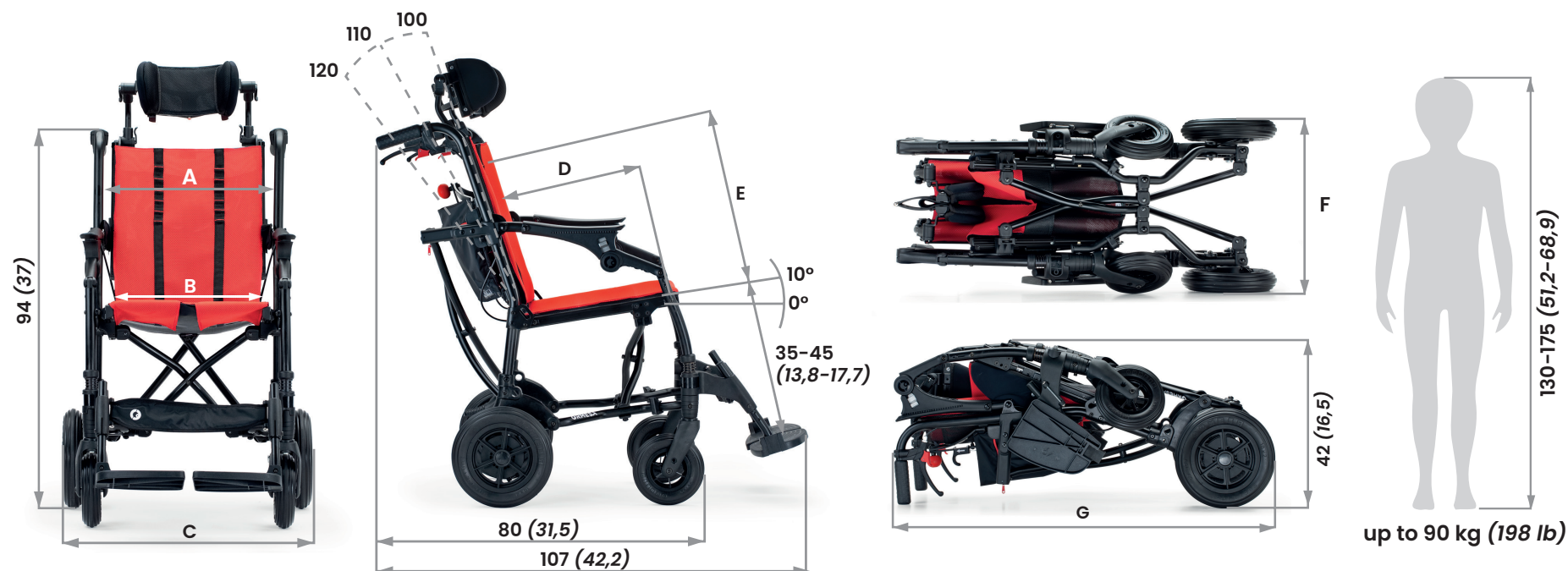

	A	B	C	D	E	F	G	Tot. weight
size 36	39 (15,4)	36 - 32 - 28 * (14,2 - 12,6 - 11 *)	59 (23,2)	36-46 (14,2-18,1)	46-56 (18,1-22)	42 (16,5)	94 (37)	15,1 (33,2)
size 40	43 (16,9)	40 - 36 - 32 * (15,7 - 14,2 - 12,6 *)	63 (24,8)	36-46 (14,2-18,1)	46-56 (18,1-22)	46 (18,1)	94 (37)	16,1 (35,4)
size 43	46 (18,1)	43 - 39 - 35 * (16,9 - 15,4 - 13,8 *)	66 (26)	36-46 (14,2-18,1)	46-56 (18,1-22)	46 (18,1)	97 (38,1)	17,1 (37,6)

* using the 958 padded pelvic side supports the width is reduced by 4 cm (1,6 in) or 8 cm (3,2 in) removing the internal padding.
36-40-43 cm (14,2-15,7-16,9 in) width without padded pelvic side supports.

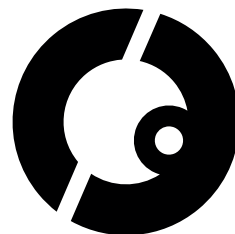

CE according to MDR (EU) 2017/745

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= **ISO 13485** =

Ormesa srl - Via A. Da Sangallo, 1 | 06034 Foligno (PG) - ITALY
Tel. +39 0742 22927 | FAX +39 0742 22637 | info@ormesa.com

www.ormesa.com