



**ORMESA®**  
• MADE IN ITALY SINCE 1980 •

**Trollino**



# Suitable for any environment







# Stability, practicality and comfort: all in one.

**Trollino** by Ormesa is a compact and light **seating system for children and teens**, ideal to move in any environment.

It is easy to handle thanks to the **international patented system** which allows you to **fold it with one single operation**. Once folded it doesn't increase its length, **occupying the minimum space**.

Trollino has a **durable and resistant structure** guaranteed by the **double cross** positioned behind the backrest and under the seat which gives greater stability and comfort.

The **upholstery is made with advanced technological materials**: it is extremely **breathable, resistant, easily washable** and **fire retardant**. It is available in **four colors**: red, blue, pink and black.

Thanks to a **wide range of available components**, it easily meets every need.



**Trollino is CRASH TESTED** according to **ISO 7176/19** and **ANSI RESNA WC/19** standards.



**TROLLINO is a MEDICAL DEVICE COMPLYING WITH THE SAFETY REQUIREMENTS DEFINED BY REGULATION (EU) 2017/745**





# Practical and light

## It fits your dynamic lifestyle



**Trollino** has an aluminium **folding frame**, painted with non-toxic epoxy powders.

**Once folded it remains compact:** this allows you **to carry Trollino** on its rear wheels like a **travel trolley**.

Its **low weight** and its **easy folding**, makes it easy to lift or move: it is perfect if you need to store it in the car or to take public transport.

The **wheels adapt to any type of terrain** offering maximum comfort to both mother and child. It is so handy that **it can be pushed with one hand**.

**Trollino** is the **ideal solution** for dynamic families often traveling and that are looking for **stability, practicality and comfort, all in one!**





# Smart upholstery

## Leading-edge cover

The upholstery of Trollino has a **honeycomb** leading-edge structure, therefore it is **extremely breathable in contact with the skin**.

The structure of the fabric creates an “air chamber” that allows both **air circulation, for maximum breathability**, and an **effective natural thermal insulation**.

**The double layer is a guarantee of:**

- › **Hygiene:** thanks to the excellent breathability, it reduces sweating and bacteria formation.
- › **Washing:** despite the upholstery is not removable, the fabric composition allows easy sanitization and washing.
- › **Quick drying.**
- › **Resistance:** High resistance to abrasion. The high resistance texture maintains its performance even at high loads (ISO 12947 test).
- › **Fire retardant:** according to EN 1021-1:2014 and EN 1021-2:2014.





# Choose your style: four diff





ifferent colors for each personality





# Technical features and adjustments



## Backrest

Reclining in three positions. Simple adjustment in height by sliding, without the risk of losing the screw. With a bag on the back to store the manual and the supplied wrench. The back cross ensures stability, comfort and endurance.



## Seat

Ergonomic and slightly tilted seat to prevent slipping of the pelvis and ensure greater comfort. Adjustable in depth by sliding and in width thanks to the additional component 958 padded pelvic side supports. Simple adjustment without the risk of losing the screw. The cross below ensures stability, comfort and endurance.





## Handles

Push handles with manual and parking brakes. Very useful handles to manage the stroller in different environments. Adjustment in height with component 951.



## Legrest

adjustable in height.

## Footrests

Tilting and swivel footrests in line with the frame without encumbrance. Adjustment in tilt with component 956.

The calf support band is also included for more comfort.



## Wheels

with "Ormesa Design" original profile.

Solid, light and **all-terrain wheels**.

**Front solid swiveling wheels:** 7" (17 cm) diameter.

**Rear solid wheels:** 10" (25 cm) diameter, with **independent drum brakes**.

**Overcoming step bar** on both rear wheels.





# Components

## The importance of an ergonomic seat

For us, **comfort and ergonomics are essential values:** Trollino makes users feel not only safe and extremely comfortable but it helps also to offer an optimal and stable posture in the most disparate circumstances.

From the backrest to the seat, passing through the headrests and the trunk side supports, **each element of Trollino is adjustable in a simple and intuitive way.** Furthermore, with its set of additional components, it offers possibilities of customization and comfort without limits.

Easy to fix and remove, handy and durable, Trollino's components are designed to the smallest detail.

The **international patented system** allows you to **fold it with one single operation, even with mounted components,** to ensure a fully satisfactory user experience.




# Components



**835 Shaped and enveloping headrest**  
adjustable in height.



**865 Padded headrest**  
adjustable in height.

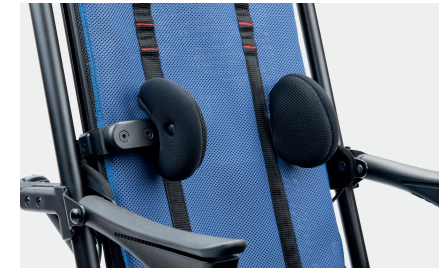


**838 Trunk side supports**  
adjustable in height, width and rotation.  
Available in two sizes.



**947 Pelvic belt with variable angle**  
Available in two sizes.




**920 Four point pelvic belt**  
Available in two sizes.

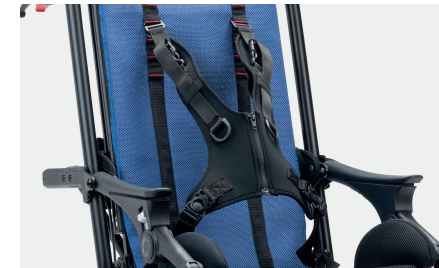

**853 Vest harness with zipper**  
Available in two sizes.



# Components

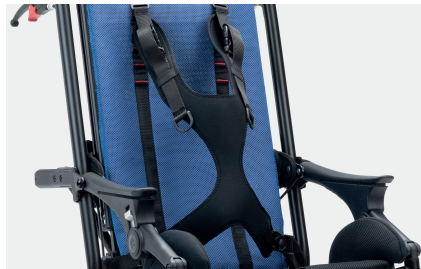

**853 slim Four point shaped harness**



**834 Padded abduction block**  
removable.



**839 Front handlebar**  
Adjustable in tilt. Easy to remove to fold the stroller.



**819 Canopy**  
Elevating, with bag and rear small window. There is no need to remove the canopy, it can be closed with the stroller!



**825 Rain cover**  
To fix on the 819 canopy.  
**Back protection Cloth (825B)**  
available for higher protection against external agents.



**818 Thermic cover**



# Components



**960 Heel rests**



**956 Footrests adjustable in tilt**



**858 Shopping basket**



**812 Front wheel directional locks**



**844 Padded side protection covers**



**951 Height adjustable push handles**

It increases the total height of max. 7 cm/2,7 in (98-105 cm/38,5-41,3 in).



**891 Set of 4 tie down hooks**


for the forward-facing transport of the pushchair in a motor vehicle (private/cars, buses..).

*See User and Maintenance Manual provided with the aid.*





# Trollino sizes and weights

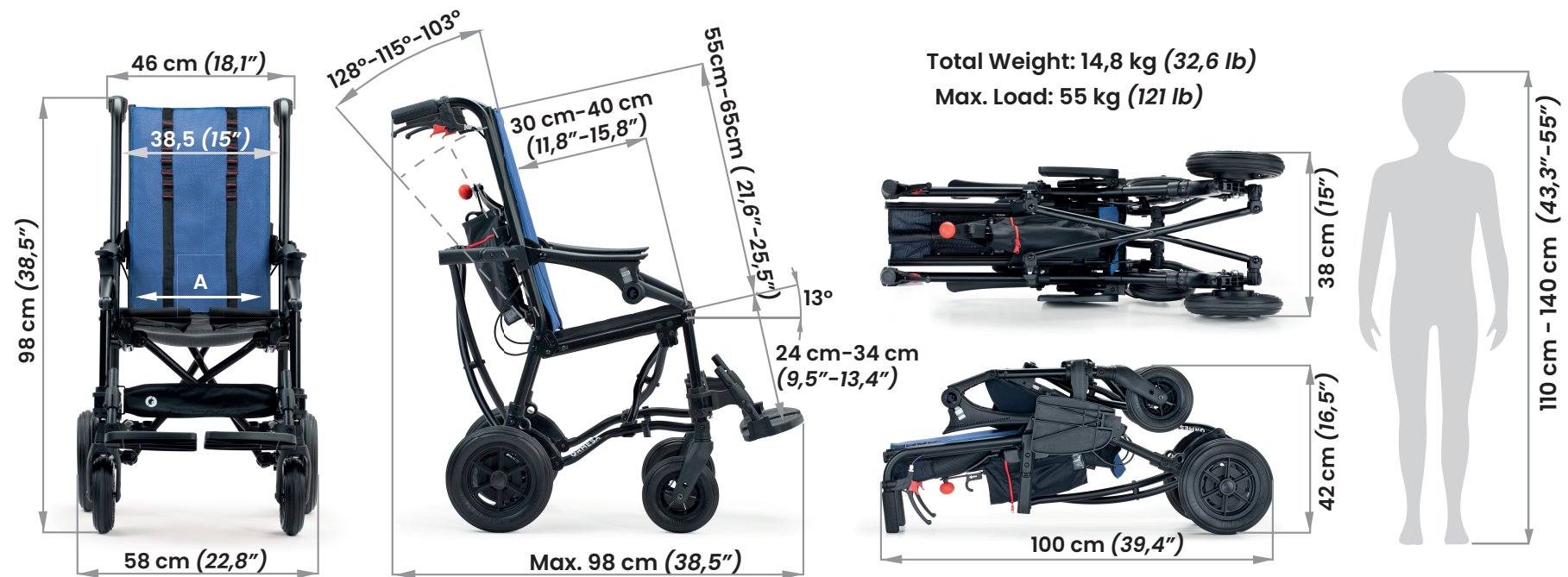

**A:** 36 cm (14"). For different width (32 and 28 cm - 12,6" and 11") please select 958 padded pelvic side supports.



Gaia and her story await  
you in the Ormesa official  
website, together with other  
protagonists of our fantastic  
global community.







**A special thanks to the families who participated with enthusiasm in our project. It is thanks to our customers and their stories that we find the energy and stimuli to evolve day by day.**

**CE** complies with MDR (EU) 2017/745

COMPANY WITH  
QUALITY SYSTEM  
CERTIFIED BY DNV GL  
= **ISO 13485** =

---

**Ormesa srl** - Via A. Da Sangallo, 1 | 06034 Foligno (PG) - ITALY  
Tel. +39 0742 22927 | FAX +39 0742 22637 | [info@ormesa.com](mailto:info@ormesa.com)

**[www.ormesa.com](http://www.ormesa.com)**