



• MADE IN ITALY SINCE 1980 •





Life made easier.

CLIP by Ormesa is a practical and fast seating system, for people with disabilities, children and adults. It is easy to handle, fold and push.

Its wheels are suitable for all terrain and **folding operations have never been so user-friendly**: CLIP is ideal for families always on the move.

Thanks to a **wide range of available components** it adapts to any needs.

Clip has a **durable structure** thanks to the **double cruise frame** positioned on the backrest and under the seat which gives more stiffness without sacrificing comfort.





CLIP is a MEDICAL DEVICE that COMPLIES WITH THE SAFETY REQUIREMENTS defined with the REGULATION (EU) 2017/745



Easy to handle

A versatile urban stroller

CLIP is available in **four sizes:** Clip 1-2-3 for children and Clip 4 for kids.

Its **low weight** and **easy folding**, make it the ideal solution for those families who are often traveling or living in the city, with obligatory and often inaccessible passages on sidewalks and crowded areas.

Padded and breathable upholstery, reinforced with internal bars, easily removable for washing.





Technical features and adjustments





Folding frame

umbrella-like closing system. It is painted with non-toxic epoxy paint.





Backrest

semirigid with inner aluminium supports: they are shapeable according to postural needs, removable to wash the upholstery. They assure the right stiffness to the backrest so that it doesn't tip over under the weight of the trunk and it also keeps the seat breathability. Three backrest positions thanks to zippers. Net basket on the back of the backrest.





Seat

inclined, breathable, semirigid seat with inner aluminium supports: the bars give the right stiffness and they help to the tightness of the seat without breaking by its use. The bars do not prevent the breathability of the upholstery.





Footrest

adjustable in height. Detachable for sizes 1-2-3 and also tilting for size 4.





Wheels

puncture-proof, suitable for all terrain.

Clip 1-2-3: front 15cm, swivelling; rear 20 cm, sliding on ball bearings with independent brakes.

Clip 4: front 17.5cm, swivelling; rear 25 cm, sliding on ball bearings with **independent drum brakes.**

Components



835 Headrest



816 Lateral supports



853 Vest harness



906 Five point harness



903 Five point vest harness



817 Leg straps for abduction

Components



894 45° pelvic belt



827 Foot straps



819 Canopy



818 Thermal cover



825 Rain cover to fix to the 819 canopy.



812 Directional locks for front wheels



Clip 1-2-3-4 sizes and weight

A: backrest width

B: seat bottom width

C: overall width

D: overall height

E: backrest inclination

F: seat bottom depth

G: backrest height

H: seat bottom inclination

I: ground to seat bottom height

L: seat to footrest height

M: handle to wheel distance

N: overall length

O: folded stroller length

P: folded stroller height

Q: folded stroller width

R: user height

Total weight
Max capacity

*seat width is 45 cm (17,7 in) if you remove the lateral paddings

clip 1

40 cm ~ 15,7 in 33 cm ~ 12,9 in 51 cm ~ 20 in 99 cm ~ 38,9 in from 0° to 12° 25 cm ~ 9,8 in 55 cm ~ 21,6 in

46 cm ~ 18,1 in from 12 to 20 cm from 4,7 to 7,8 in

73 cm ~ *28,7 in*

80 cm ~ 31,4 in 110 cm ~ 43,3 in

30 cm ~ 11,8 in

34 cm ~ *13,3 in*

from 110 to 125 cm from 43,3 to 49,2 in

10,8 kg ~ *23,8 lb*

30 kg ~ 66,1 lb

clip 2

40 cm ~ 15.7 in 33 cm ~ 12.9 in 51 cm ~ 20 in 99 cm ~ 38,9 in from 0° to 12° 30 cm ~ 11,8 in 59 cm ~ 23,2 in 180 48 cm ~ 18.8 in from 21 to 31 cm from 8.2 to 12.2 in 73 cm ~ 28.7 in 88 cm ~ 34.6 in 110 cm ~ 43,3 in 30 cm ~ 11,8 in 34 cm ~ 13.3 in from 120 to 140 cm from 47,2 to 55,1 in 11 kg ~ 24,2 lb

clip 3


40 cm ~ 15.7 in 36 cm ~ 14,1 in 51 cm ~ 20 in 104 cm ~ 40,9 in from 0° to 12° 35 cm ~ 13.7 in 64 cm ~ 25,1 in 18° 49 cm ~ 19.2 in from 27 to 37 cm from 10.6 to 14.5 in 79 cm ~ 31.1 in 96 cm ~ 37.7 in 115 cm ~ 45.2 in 32 cm ~ 12,5 in 35 cm ~ 13.7 in from 135 to 155 cm from 53.1 to 61 in 11,2 kg ~ 24,6 lb 40 kg ~ 88,1 lb

clip 4

45 cm ~ 17.7 in 33 cm* ~ 12,9 in* 57 cm ~ 22,4 in 113 cm ~ 44,4 in from 0° to 10° 38 cm ~ 14.9 in 71 cm ~ 27,9 in 20° 59 cm ~ 23.2 in from 37 to 47 cm from 14,5 to 18,5 in 82 cm ~ 32.2 in 105 cm ~ 41.3 in 115 cm ~ 45.2 in 56 cm ~ 22 in 34 cm ~ 13.3 in from 150 to 175 cm from 59 to 68,8 in 13,5 kg ~ 29,7 lb 76 kg ~ 167,5 lb



35 kg ~77,1 lb

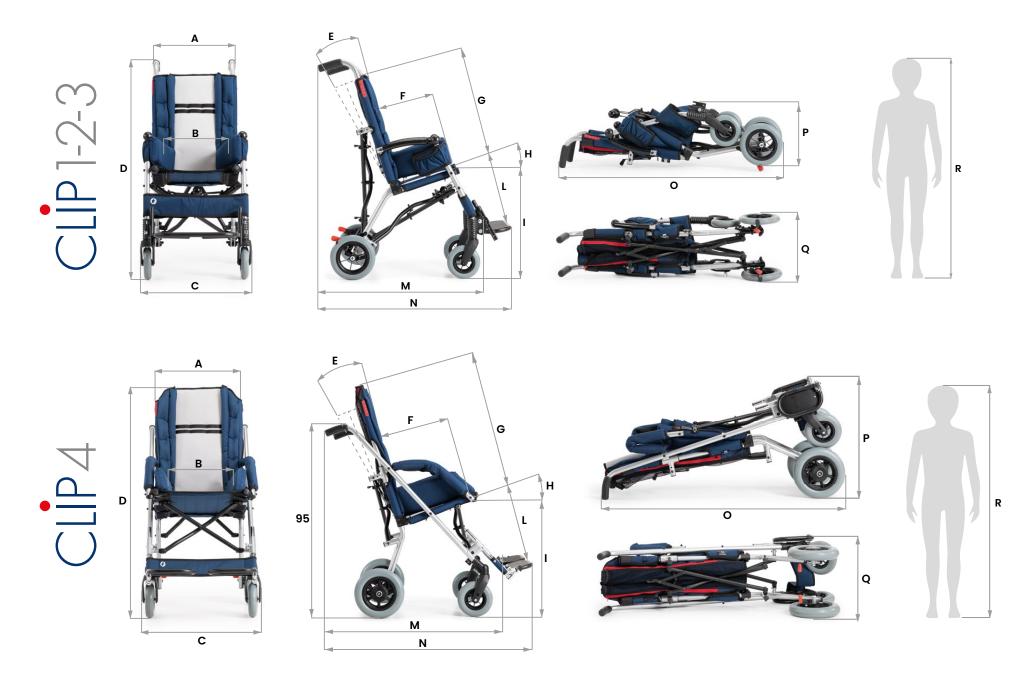



(€ complies with MDR (EU) 2017/745

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 13485 =

Ormesa srl - Via delle Industrie 6/8 | 06034 Foligno (PG) - ITALY Tel. +39 0742 22927 | FAX +39 0742 22637 | info@ormesa.com

www.ormesa.com