

# Flexmobil i6 Comfort II

With the introduction of Flexmobil Comfort II, we are further enhancing the comfort of people who are sitting in their electric wheelchairs more frequently and for longer periods of time. Flexmobil Comfort II is an electric wheelchair for indoor use, affording users a new level of freedom of movement in the home and at work. The chair features an innovative chassis with a more generous seat height and lower noise level. The new Comfort II is an evolution of the previous Comfort wheelchair, with a new, sliding back support function that delivers an enhanced level of seating comfort thanks to the introduction of new cushion technology. The chair is also easier to adjust.

Type of impairment


minor  moderate  extensive

indoor  outdoor  child  adult



## SEAT SYSTEM

Comfort

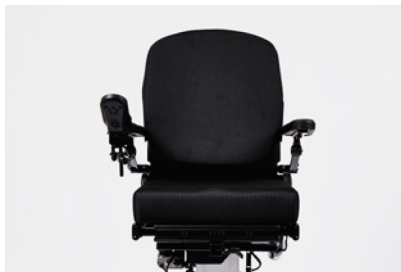

## UPHOLSTERY



Back fabric

The comfortable Comfort seat provides ergonomic support for the entire body, allowing you to sit comfortably the whole day. The seat is available in three widths, with a low or high back support. You can also choose either the Base model, which is built up in multiple foam rubber layers, or the Contour, which combines multiple foam rubber layers with other elements of memory foam. The position of the seat is fully adjustable up to 60 cm high, allowing users with longer legs to find a comfortable sitting position.

Seat (width x depth)	40x40, 40x45, 40x50, 40x55, 45x40, 45x45, 45x50, 45x55, 50x40, 50x45, 50x50, 50x55 cm
Seat width (between armrests)	33-73 cm
Seat depth	30-60 cm
Seat height (height to underneath the seat)	39-64 cm, 48-88 cm
Seat tilt, electrical, crank & turnbuckle	-12° till 20°
Back (width x height)	40x45, 45x45, 50x45, 40x54, 45x54, 50x54 cm
Back tilt	22°
Upholstery	Black fabric
Armrest, height, adjustable	10-39 cm
Armrest platform	30x8x3 cm

The seat height is adjusted electrically and is available in high and low versions. Both are listed after each other in the technical data in the order "low, high".

## CHASSIS

Flexmobil i6



Flexmobil i6 makes it easy to both "sit near the floor" or "stand up and work" thanks to the generous min and max positions. The chassis is available with two different actuators. Thanks to the central placement of the drive wheel, the Flexmobil i6 is very easy to drive. The chair turns on its own axis, with an extremely small turning circle that doesn't need much space.

Chassis, width	57 cm
Chassis, length	80 cm
Total weight	110 kg
User weight, max.	150 kg
Turning circle	55 cm
Max speed	4 km/h
Max operation time	8,5 h
Electronics	Dynamic/PG
Obstruction clearance	3 cm
Drive wheel size	225 mm
Link wheel size	125 mm